

PRESORT STANDARD
U. S. POSTAGE
PAID
PERMIT NO. 3280
SALT LAKE CITY, UT

UTAH MEDICAL PRODUCTS INC.
7043 South 300 West • Midvale, Utah 84047
800-533-4984 • 801-566-1200
www.utahmed.com



NOT ALL
NEONATAL PICCS
ARE CREATED
EQUAL ...

PICC-NATE[®] : CREATED TO BE DEVELOPMENTALLY FRIENDLY IN MEETING YOUR FRAILEST PATIENTS' NEEDS

Utah Medical Products, Inc. (UTMD) Gesco neonatal products have a proven reputation for high quality of construction. UTMD designed PICC-Nate exclusively for the unique care criteria of neonates:


- ◆ Barium sulfate-loaded catheter provides radiopacity for maximum visualization while maintaining the safety and effectiveness of the catheter.
- ◆ Silicone catheter is non-reactive to body tissues and fluids, non-supportive of bacterial growth, non-thrombogenic, and minimally irritating to vessel walls and organs.
- ◆ The inherent properties of silicone allow the catheter to remain in the center of the vein minimizing damage and inflammation to the vessel walls.

UTMD's neonatal products should be integral to any NICU's Developmental Care program.

The low-profile hub reduces tissue injury. Less trauma to a baby's skin can reduce the number of necessary dressing changes, and allow the PICC to remain indwelling for a longer duration.

The extended hub option eliminates the need for an extension set and distances the bulkiest part of the catheter from the delicate skin of the baby, thereby minimizing skin irritation.

The suture wing design and tape notches are profiled for ideal hub securement with minimal skin irritation and easy site maintenance.

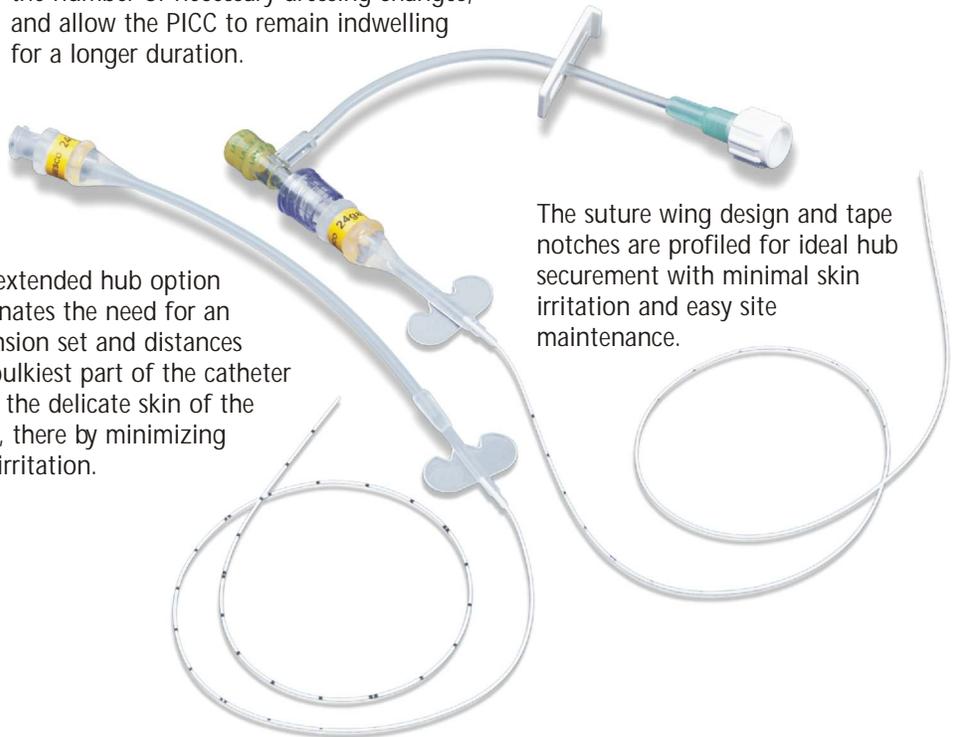

Funneled entry promotes fast PICC introduction.

"Quick Flashback" cannula indicates immediate verification of vessel penetration.



SILICONE	1.9 Fr (24 Ga)	3.0 Fr (20 Ga)
Catheter & Introducer	PB-2PS	PB-3PS
Complete Tray	PC-2PS	PC-3PS
SILICONE w/ EXTENDED HUB	1.9 Fr (24 Ga)	3.0 Fr (20 Ga)
Catheter & Introducer	PB-2LPS	PB-3LPS
Complete Tray	PC-2LPS	PC-3LPS
TECOFLEX POLYURETHANE	1.9 Fr (24 Ga)	3.0 Fr (20 Ga)
Catheter & Introducer	PB-2PP	PB-3PP
Complete Tray	PC-2PP	PC-3PP
INTRODUCER ONLY	PEEL2 (22 ga)	PEEL3 (19 ga)

For more information or to place an order, visit:

www.utahmed.com

or call:

1-800-533-4984